

BINGO EVOLUTION

Innovative and modern rehab pushchair.



Designed and hand built in Germany.

BINGO EVOLUTION



BINGO Evolution - A good product is constantly improving.

Pictured configuration:

Mobility base swivelling front wheels with polyurethane wheels (flat-free) on black rims. Seat unit with upholstery in "turquoise soft" including laser cut side covers matching the upholstery colour. Headrest pads (low version). Hip supports and thigh supports. Frame suspension system, spring set matching the upholstery colour. Foot plate in boat form. Soft accessory bag. Frame colour: "silver".



Innovative and modern rehab pushchair.

Rehab pushchairs are an essential part in children-centred care. With individual solutions and modern designs **HOGGI** is setting standards since 2005.

HOGGI philosophy. Passion and continuous product optimising.

The result of our own ambition of a continuous further development and product optimising has a name: **BINGO Evolution.**

A good product is constantly improving.

Profit by our longtime experience and our passion.

Cool design. HOGGI laser cut side covers.



HOGGI laser cut side covers.



Compact. Light. Sporty.

Modern, compact and sporty design.

BINGO Evolution - an impressive design with modern, sporty and compact appearance. Select your favourite combination from five upholstery colours and two frame colours and create your own **BINGO Evolution**.



Multifunctional seat unit. Adjustment possibilities.

45° tilt in space.

The new construction enables a quick and easy tilt in space of the seat unit up to 45° in 4° steps.

Recline.

The recline is also easily and quickly variable in 7,5° steps.

High seat comfort thanks to tilt in space and recline.





Central positioned. Easily operated.

Operating elements of tilt in space and recline are positioned centrally on the backrest and offer an easy and safe "handling".



Quick rotation of the seat unit. Without tools.

According to the individual need of the child, eye contact can be kept or the seat unit is oriented in driving direction. The weather situation can also necessitate a quick rotation of the seat unit.

You can still profit of the function of tilt in space and recline even if the seat unit is inserted opposite to the driving direction.

Frame suspension system.

The frame suspension system which is four times variabel in a comfortable way absorbs shocks from the ground and reduces the risc of spasms. We offer two suspension systems in two intensities.

Rotation of seat unit. Frame suspension system.




Seat unit in driving direction.




Release button.



Seat unit against driving direction.

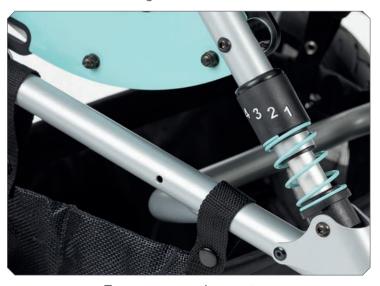

Frame suspension system.

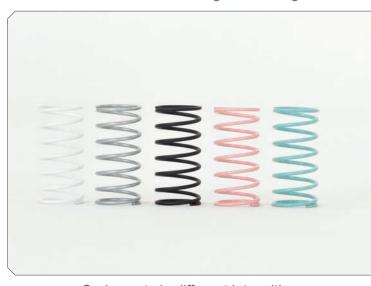

Spring sets in different intensities.



Comfortable release of folding mechanism on push bar.

The release buttons for the innovative folding mechanism are positioned on the push bar, comfortable to reach and easy to release. The seat unit folds as well.

Extreme compact folding size.

For transportation or car ride only minimal manipulations are necessary to fold and lock your **BINGO Evolution.** The buildup is just as easy.

Innovative folding mechanism. Extreme compact folding size.



Height adjustable push bar.



Compact folding size.



Variations of foot plate.



Wheel lock.



Headrest pads. Trunk supports. Thigh supports.

The headrest pads (low version), the trunk supports as well as the thigh supports ensure a sporty and filigree optic. *HOGGI* offers a wide range of supports for perfect patient-centred care. All supports are variable in width, depth and angle. Please feel free to ask us!

IUS. Individual upholstery system. Different upholstery colours.

Customisable upholstery including *HOGGI* laser cut side covers matching the upholstery colour. Select your favorite colour from five different upholstery colours.

Supports. Upholstery. Side covers.





IUS. Individual upholstery system.

- 1. Back cushion
- 2. Cover for trunk supports
- 3. Seat cushion
- 4. Side covers


BINGO EVOLUTION

Growth capacity.

BINGO Evolution has a high growth capacity in back height, seat depht, seat width and lower leg lenght.

Individual adjustability.

The multifunctional seat unit offers individual adjustment possibilities for your child thanks to a wide range of supports.

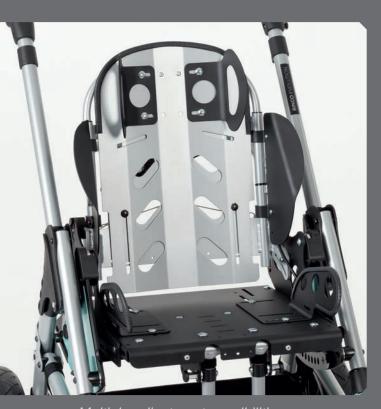


Multiple adjustment possibilities.



Huge growth capacity.

Central operating elements.

- 1. Release for tilt in space
- 2. Release for recline

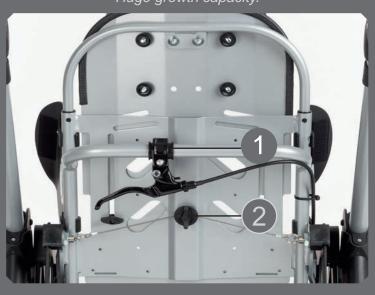

BINGO MOBILITY Base

BINGO Mobility base is specially designed for the use of seat shells:

- Compact folding size
- Tilt in space 45°
- Quickly change the seat from front to rear facing



BINGO Mobility frame.



Adjustment range tilt in space.



Adjustment range recline.



Seat frame & back base.





BINGS EVOLUTION

Indoor bases COBRA and BASIC.

COBRA has been specially designed as an indoor base for the **BINGO Evolution** seat unit with 45° tilt in space and continuously height adjustment by gas pressure springs for various indoor applications:

- as indoor base for the BINGO Evolution seat unit
- as mobility base for orthopedic seating systems and seat shells

BASIC is a classic indoor base for the **BINGO** Evolution seat unit.



COBRA.



With seat frame & back base.



COBRA, lowest position.



BASIC.



Height adjustable range.



With seat frame & back base.





Specifications	Size 1	Size 2	Size 2 XL
Seat width	18 - 30 cm	23 - 35 cm	23 - 35 cm
Seat depth	16 - 34 cm	22 - 40 cm	30 - 45 cm
Tilt in space seat unit In driving direction Opposite to the driving direction	45° 45°	45° 45°	45°
Back height	45 - 63 cm	54 - 75 cm	62 - 83 cm
Recline seat unit Recline seat unit Recline seat unit	90° - 160° 90° - 160°	90° - 160° 90° - 160°	90° - 135° -
Lower leg length	16 - 37 cm	16 - 40 cm	16 - 40 cm
Height of push bar	75 - 110 cm	75 - 110 cm	86 - 118 cm
Overall width	63 cm	68 cm	68 cm
Folding size (L x W x H)	78 x 62 x 38 cm	78 x 67 x 38 cm	78 x 67 x 38 cm
Weight limit Seat unit Mobility base	50 kg 60 kg	50 kg 60 kg	50 kg 60 kg

Safty aspect: Successfully crashtested.

BINGO Evolution is crashtested successfully according to ISO 7176-19. **BINGO Evolution** is admitted for use as a seat in a motor vehicle.



HOGGI

Partnering together to mobilise kids.

HOGGI GmbH

♥ Eulerstraße 27 D-56235 Ransbach-Baumbach

Fon: (+49) 2623 / 92 499-0 Fax: (+49) 2623 / 92 499-99

www.hoggi.de