

www.hoggi.de

Made in Germany
by HOGGI 

ZIP

Solid and compact rehab buggy.



Designed and hand built in Germany.

ZIP

Solid and compact rehab buggy.



ZIP®

Pictured configuration:

Solid and compact rehab buggy made of aluminium. Frame colour: „black“. Reinforced black nylon cover. Footrest: infinitely variable in height, angle adjustable and swingaway. Wheel lock. Bearing wheels with polyurethane tires on black rims. Swivel lock.

ZIP

Modern. Sporty. Extreme compact. Light.

Solid, compact and light rehab buggy made of aluminium:


- Modern and sporty design
- Easy handling and extreme compact folding size
- Innovative accessories for perfect, individual patient-centred care
- High usability due to huge growth capacity
- Available in different colours for frame, nylon and upholstery



Modern.



Sporty.

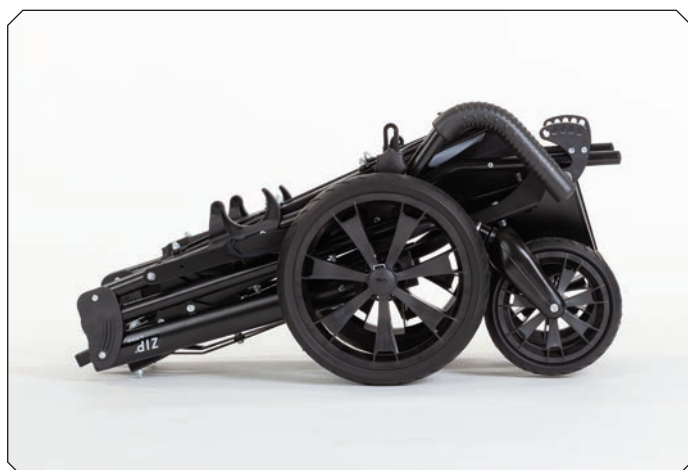

Extreme Compact.

ZIP

ZIP® your buggy. Always flexible - in any weather.

! PU chrome rims are no longer available.



Turquoise Soft.



Berry Soft.



Grey Soft.



With winter warmer.



With foldable canopy.



With canopy and rain cover.

Handy in workaday life.

Parents and carer have specific requirements on rehab buggies for their children. They require a reliable product giving excellent mobility and independence - the ZIP® rehab buggy will meet all these requirements. With its unique compact folding size of only 78 cm (91 cm ZIP® size 2) the ZIP® rehab buggy fits in most car boots and lives up to its name.

ZIP

Solid and compact rehab buggy.



„COLOUR-LINE Berry Soft“.

Pictured configuration:

Solid and compact rehab buggy made of aluminium. Frame colour: „black“. Reinforced black nylon cover. Foldable backrest and seat upholstery. Upholstery colour: „Berry Soft“. Footrest: infinitely variable in height, angle adjustable, foldable and swing away. Wheel lock. Bearing wheels with polyurethane tires on black rims.

ZIP

Infinitely variable footrest. Optional accessories.



Aluminium footrest.



Swivel lock.

Perfect patient-centred care. Chic and elegant solutions.



Grab rail.



4-point lap belt.



Lap belt.

Useful and comfortable in basic equipment. Innovative optional accessories.

The footrest infinitely variable in height and angle adjustable as well as the attendant brake are included in the ZIP® basic model. Various sophisticated accessories (22!) give the opportunity to customise your buggy according to patient-centred care. ZIP® keeps its straight lines even „fully equipped“!

All ZIP® accessories can be found with a foto in the interactive order form or online in the ZIP® accessories gallery at : www.hoggi.de.



Solid and compact rehab buggy.

Specifications	Size 1	Size 2
Seat width	35 cm	40 cm
Seat depth	30 or 35 cm	35 or 40 cm
Seat angle	18°	18°
Height push handle	110 - 112 cm	110 - 112 cm
Backrest height		
Backrest, standard	60 cm	68 cm
Backrest, long	66 cm	
Backrest angle	90° - 100°	90° - 100°
Lower leg length	20 - 39 cm	25 - 44 cm
Footrest angle	-10° to + 10°	-10° to + 10°
Wheel diameter		
Front wheel	20 cm	20 cm
Rear wheel	30 cm	30 cm
Overall width	59 cm	65 cm
Folding size	78 x 36 x 39 cm	91 x 44 x 40 cm
Load capacity	50 kg	60 kg
Weight	12,5 kg	13,5 kg
Turning radius	126 cm	140 cm

Safety aspect: Successfully chrashtested.

ZIP® is chrashtested successfully according to ISO 7176-19.

ZIP® is admitted for use as a seat in a motor vehicle.

The belt connection points for the transport in a motor vehicle are included in the basic model!





Solid and compact rehab buggy.

HOGGI

Partnering together to mobilise kids.

HOGGI GmbH

 Eulerstraße 27
D-56235 Ransbach-Baumbach

 Fon: (+49) 2623 / 92 499-0
Fax: (+49) 2623 / 92 499-99

 info@hoggi.de

www.hoggi.de