

www.hoggi.de

Made in Germany
by HOGGI 

CESA

Sporty active-wheelchair for kids & adolescents.



Designed and hand built in Germany.

CESA

Product information.

CESA - sporty wheelchair for kids & adolescents.

With **CESA HOGGI** launches a new sporty wheelchair and extends the range of active wheelchairs for kids and adolescents. The **box frame construction** provides **CESA** with a huge rigidity despite its low weight. Thanks to the special box frame construction seat height, seat angle, wheel base and active degree can be adjusted independently of the chassis.

An expensive and defective calibration/adjustment of wheel fork and wheel track of the rear wheels is prevented completely.

The **patented frame concept** consists of sparse components and can be adjusted only by one hexagon allen wrench.

Further USP:

The **wheel base** can be adjusted independently of the active degree with the aid of a **slide axis**. In this way the wheel base can be adapted to the seat depth and it is possible to keep a steady footrest angle sovereign of the seat depth.

In addition the weight can be scattered optimally on rear and front wheels - an important feature for energy efficiency and dynamic driving quality.


CESA convinces also with an enormous growth capacity:

- Seat depth up to 10 cm without any application of additional items
- Seat width up to 2 cm with distance elements

In addition to many colour and design variants, **CESA** offers a wide range of options and accessories to meet both individual design requirements and sophisticated care concepts.



Box frame construction.

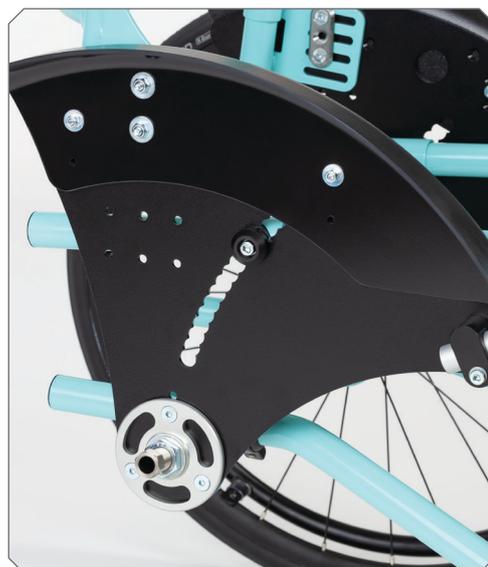

Side plate (seat angle adjustment).



Modern, sporty design.

CESA

Specifications.

Specifications	Frame size 1	Frame size 2	Frame size 3
Seat width	24 - 30 cm	28 - 36 cm	34 - 44 cm
Seat depth	24 - 36 cm	26 - 42 cm	32 - 48 cm
Seat height	36 - 42 cm	37 - 45 cm	41 - 49 cm
Seat angle (adjustable)	0° - 12°	0° - 12°	0° - 12°
Backrest height	20, 25, 30, 35 cm	20, 25, 30, 35 cm	20, 25, 30, 35 cm
Backrest angle			
Rigid back	-10° up to +10°	-10° up to +10°	-10° up to +10°
Reclinable back	80° - 105°	80° - 105°	80° - 105°
Lower leg length	10 - 39 cm	10 - 43 cm	10 - 47 cm
Footrest angle	-10° up to + 10°	-10° up to + 10°	-10° up to + 10°
Rear wheel diameter	20"	22"	24"
Front wheel diameter	4", 5", 5,5"	4", 5", 5,5"	4", 5", 5,5"
Camber	3°, 6°, 9°, 12°	3°, 6°, 9°, 12°	3°, 6°, 9°, 12°
Weight limit	60 kg	80 kg	100 kg



Back base, rigid, angle adjustable.



Back base, angle adjustable, foldable.



CESA with seat panel.



Sporty active-wheelchair for kids & adolescents.

HOGGI

Partnering together to mobilise kids.

HOGGI GmbH

📍 Eulerstraße 27
D-56235 Ransbach-Baumbach

☎ Fon: (+49) 2623 / 92 499-0
Fax: (+49) 2623 / 92 499-99

✉ info@hoggi.de

www.hoggi.de